*No. 1 choice  
of hospitals  
& mothers*

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## **Instructions for Use**

	<b>Page</b>
Various helpful information	3 – 5
Cleaning/Sterilization	5 – 6
Instructions for use	7 – 10
Troubleshooting	11
Parts list	12
Cooler Bag	13
Technical Data	14

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## **Breastmilk is best for your baby**

The close bond between mother and child is promoted in a unique way by breastfeeding. Breastmilk is the only source of nutrition that contains all the necessary elements such as protein, fat, sugar, etc. in the correct quantities and composition. The antibodies contained in breastmilk are like an injection for your child against infections, diarrhea, colic and otitis media.

The benefits for the mother include contraction of the uterus in a natural way, reduced risk of breast and ovarian cancer and, in most cases, weight reduction. If breastfeeding is not always possible, a Medela breastpump can ensure that your child gets the valuable breastmilk it needs. Our best wishes to you and your baby.

## **Hygiene**

We recommend that the hands be thoroughly washed with soap and water before using the breastpump. The breast and surrounding area should be cleaned with a damp cloth. Use of soap or alcohol can dry out the nipple skin and cause cracked nipples.

If the collected breastmilk is used to feed a premature or sick baby, it is important that the bacteria count is kept as low as possible. This is accomplished by stripping the breast: hand express the first few squirts of breastmilk and discard the milk.

# General Information

Be comfortably seated and relax. Use a footrest, when pumping your breast or when feeding your baby. Elevating your feet puts you into proper sitting position to avoid stress on the lower back and perineum. Hold the breastshield between index finger and middle finger and press it gently but firmly over the nipple, while supporting the breast a little with the same hand.

Experiment with different vacuum levels. Choose the level which works best for you and does not cause you any discomfort. Some time may be needed before milk flow is adequate. This is normal and the quantity of milk produced may vary.


## Pumping Time

Milk is produced on a supply and demand basis. Frequent pumping and completely emptying the breasts encourages the increase of milk supply. The milk ejection response can be stimulated through breast massage before or while pumping. Increasing the frequency of pumping sessions is more effective for increasing the volume of milk than increasing the duration of each episode. If the breasts are not engorged and if not otherwise directed by the doctor, use the table below as a guide.

### Single Pumping:

#### First breast

5 – 7 minutes  
3 – 5 minutes  
2 – 3 minutes
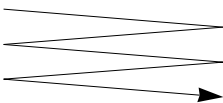
#### Second breast

5 – 7 minutes  
3 – 5 minutes  
2 – 3 minutes

**Double Pumping:** Pumping time will be 10 to 15 minutes.

## Storage of Breastmilk

Do not let breastmilk stand at room temperature. To inhibit the growth of bacteria, breastmilk which is not used immediately must be refrigerated or, if stored for more than 24 hours, must be frozen.

For freezing, do not fill bottles more than  $\frac{3}{4}$  full, to allow space for expansion or use the **Pump & Save™ Breastmilk Bags**. In the deep freezer (–18 °C) breastmilk can be safely stored for approximately three months. Never add warm milk to frozen milk.

Thaw frozen milk by holding the container under warm water or leave it for 12 hours in the refrigerator. Do not thaw frozen breastmilk in a microwave or in boiling water. Before feeding, shake the bottle to blend any fat that has separated. Check with your hospital for specific storage instructions.

## Cleaning/Sterilization

### **Important:**

For cleaning/sterilizing purposes, the valve must always be removed from the breastshield and the membrane is to be taken off the valve head. For hygienic reasons it is important to rinse the breastshield from all directions, including upside down. Wipe the breastpump unit with a damp (not wet) cloth. If the tubing show traces of milk, they have to be cleaned as well as the gaskets. Remove the gasket according to the description on page 11. Wash gaskets and tubing with warm soapy water, rinse with clear water and dry thoroughly.





**Please note that this material is not suitable for sterilization in a microwave oven.**

**Prior to initial use and subsequently once daily,** the individual components which come in contact with milk must be rinsed and placed in a saucepan of sufficient capacity. Fill the latter with cold water until the components are well covered, heat it up and allow to boil for 5 minutes.



After draining the water off, the hot items should be spread out on a clean towel to cool and dry. Instead of boiling, you can put the disassembled parts in a dishwasher and run a complete standard program.

**Each time they have been used,** all items (except for the breastpump) should be washed in lukewarm water containing a gentle cleaning agent and then rinsed in fresh water. Allow parts to dry in the air on a dry cloth.

## Instructions for use

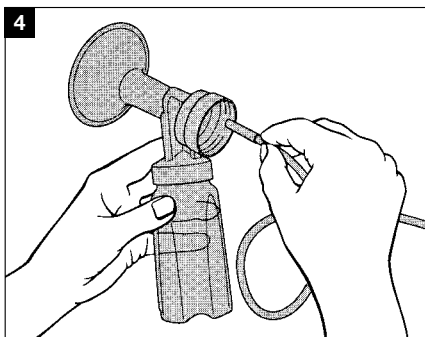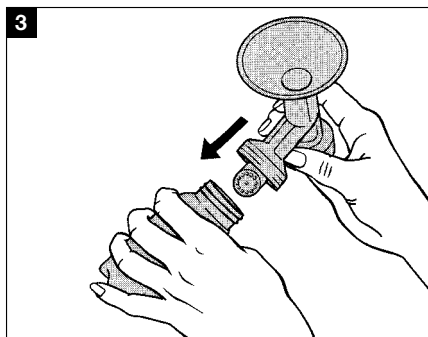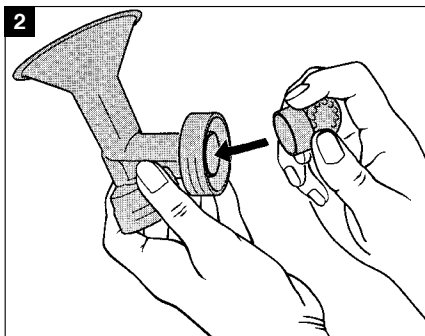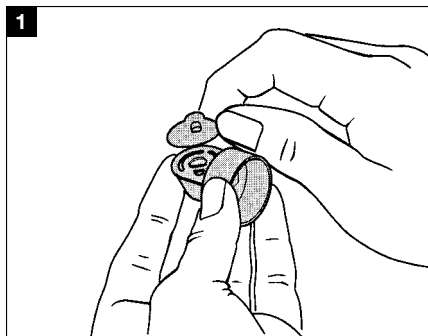
1. Snap the membrane onto the valve head.  
(Remove membrane for cleaning.)
2. Fit the valve to the breastshield (do not push too hard, so that it can easily be removed for cleaning).
3. Screw the bottle onto the breastshield.
4. Insert the end of the tubing with the white connector into the breastshield.

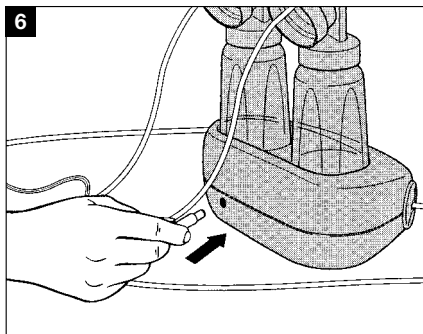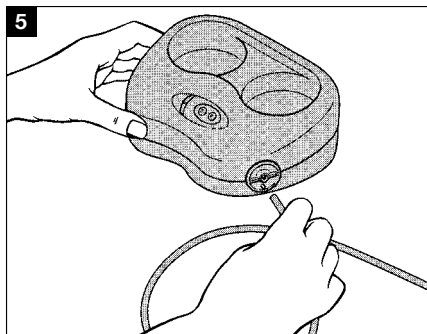5. Connect the other end of the tube to the breastpump.

Follow steps 1–5 for the second breastshield.


6. Insert the power connector into the connection socket on the breastpump. Plug in the mains adapter. **Use only the original Medela adapter. Other adapters may damage the pump or cause problems during use.**

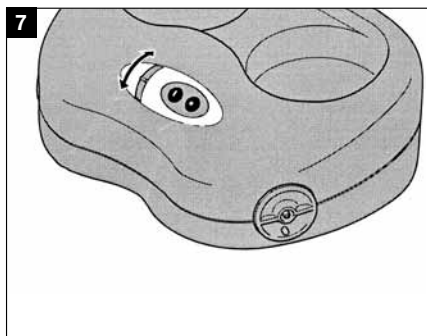7. Set the vacuum regulator to minimum and switch on the breastpump. Adjust the vacuum to your **Maximum Comfort Vacuum™**. This is the highest vacuum you can use and still feel comfortable (different for every mother). To find your Maximum Comfort Vacuum increase vacuum until pumping feels slightly uncomfortable (not painful), then decrease slightly.

**Research has shown that you will pump more efficiently – meaning get more milk in less time – when pumping at Maximum Comfort Vacuum™.**

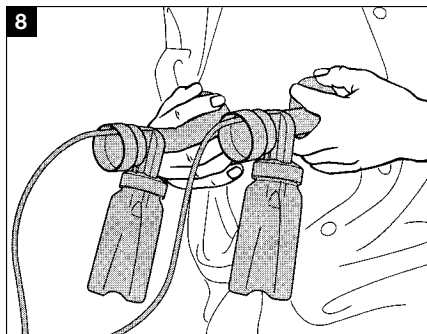8. Center the nipple in the breastshield tunnel. Turn on the pump and check the following points: Your breastshield is the right size if your answer to the following questions is «Yes»:



- Is your nipple moving freely in the tunnel?
- Is minimal or no areolar tissue being pulled into the tunnel of the breastshield?
- Do you see a gentle, rhythmical motion in the breast with each cycle of the pump?
- Do you feel the breast emptying all over?
- Is your nipple pain-free?

If your answer to any of the above questions is «No», try a larger (or smaller) breastshield.


A correctly fitted breastshield avoids compression of any milk ducts during pumping to ensure optimal breast emptying and maximum milk output.

## Single pumping

Follow steps 1 to 8, pages 7–9.

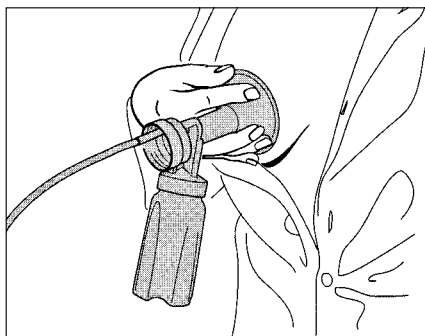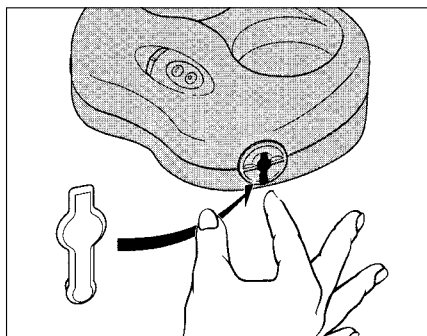Close the second connection on the breastpump with the closing cap.
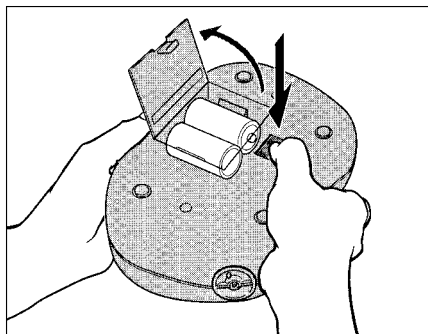
## Battery operation

If the mains adapter cannot be used, the Mini Electric Plus Breastpump can also be used with batteries. Two C Baby LR 14 1.5 V alkaline batteries are required for this purpose. Leak-proof batteries should always be used to protect the unit. If the suction rhythm becomes too slow, replace the batteries.

- Press the tab (PRESS) to open the cover of the battery compartment.
- Insert two 1,5 V batteries (c.f. «Battery operation» for the required battery specifications). **Attention:** Observe correct polarity (+ and -) as indicated inside the battery compartment.
- Close the battery compartment until the cover snaps into place.



The batteries do not need to be removed if the mains adapter is used for pumping.

## Troubleshooting


**If you do not feel any vacuum at the breast or if the pumping action is too weak:**

- Check that the breastshield forms a complete seal around the breast.
- Check that the valve is clean and properly fitted (steps 1 and 2, page 7). If necessary, remove it for cleaning/sterilization and reassemble it.
- Check that the tubing are connected correctly.
- For single pumping, check whether the second connection on the pump has been closed.

## Changing the gaskets

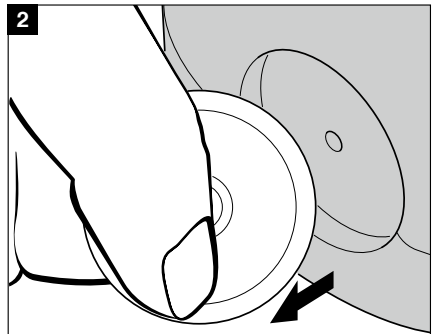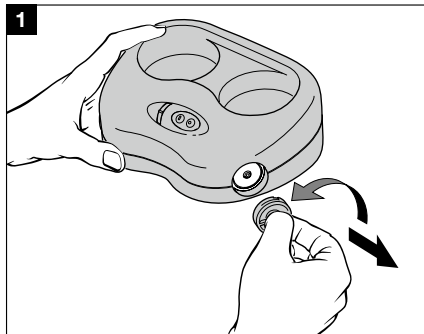1. Turn the cover counterclockwise and pull it out.
2. Remove the gasket for cleaning (see page 5).
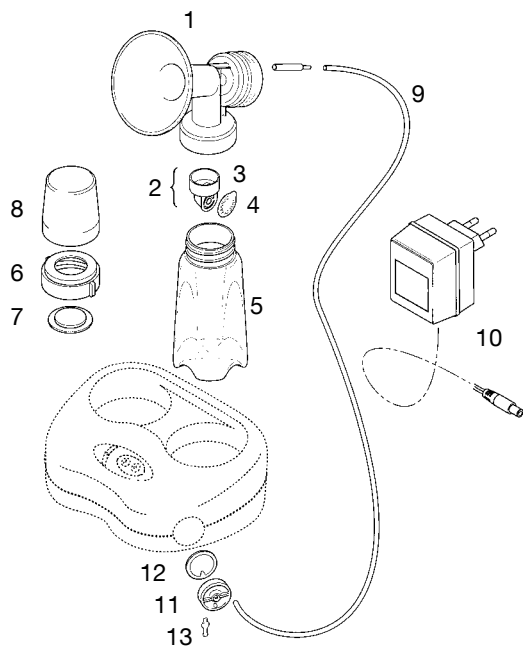
## Parts List

1. Breastshield
2. Valve (2-part)
3. Valve head
4. Valve membrane
5. Bottle
6. Lid
7. Disk
8. Cap
9. Silicone tubing with connector
10. Mains adapter for 230 V
11. Cover
12. Gasket
13. Plug to cover
14. Instructions for use

## Order No.

- 810.0131
- 800.0423
- 810.0475
- 820.7013
- 811.0020
- 810.7062
- 810.0473
- 810.0474
- 800.0436
- 920.0015
- 810.0447
- 820.0084
- 810.0332
- 190.0714

*Does not contain Bisphenol-A, DEHP (Phthalate) or Latex.*

*Subject to change.*

*Content may vary from picture.*

## **CityStyle™ Breastpump Bag**

A stylish and attractive bag to transport a Medela breastpump and expressed breastmilk. Offers additional room for the breastpumps Swing, Harmony, Mini Electric Plus, Mini Electric or Manual Breastpump.

Order-Nr. 200.0625

## Technical Data

Suction strength:	Infinitely-variable up to approx. 220 mmHg
Suction cycles:	approx. 34 cycles/minute
Suction unit:	3 VDC/500 mA or 2 x 1,5 V batteries (Two C Baby LR 14, 1,5 V alkaline)
Mains Adapter:	Input: 230 V/50 Hz, Output: 3 VDC/500 mA

### Environmental Conditions for Transport and Storage:

– Temperature:	–20 °C to +50 °C
– Relative Humidity:	20% to 95%
– Atmospheric Pressure:	70 kPa to 106 kPa

Symbols:

Protection Type B

Protection Class II (double insulation)

I / O

ON / OFF

## Warranty

1-year warranty for the unit from the date of purchase, excluding pump kit. The manufacturer is not liable for any damage or consequential damage caused by incorrect operation, inappropriate usage or use by unauthorized persons.

## Disposal

The device should not be disposed of as municipal solid waste. Use your local collection and return systems for old electronic equipment. If old electronic equipment is not disposed of correctly, this can have harmful effects on the environment and on public health.

CE-mark according to 93/42/EEC

0123

## Note

Use the breastpump in dustfree surroundings only. Unit not to be left or to be used near solar radiation or other heating unit. Equipment not suitable for use in the presence of a flammable anaesthetic mixture. Separation from mains power supply is achieved by disconnecting the mains adapter.

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Medela UK Ltd.